



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## 5-Methoxysalicylic acid

Cat. No.:	HY-W007856		
CAS No.:	2612-02-4		
Molecular Formula:	C <sub>8</sub> H <sub>8</sub> O <sub>4</sub>		
Molecular Weight:	168.15		
Target:	Endogenous Metabolite		
Pathway:	Metabolic Enzyme/Protease		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	6 months
		-20°C	1 month

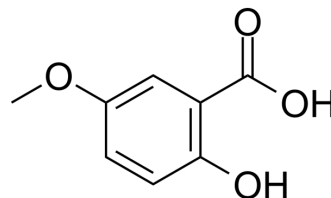

### SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (594.71 mM; Need ultrasonic)			
		Solvent Concentration	Mass	
			1 mg	5 mg
	Preparing Stock Solutions		10 mg	
	1 mM	5.9471 mL	29.7354 mL	59.4707 mL
	5 mM	1.1894 mL	5.9471 mL	11.8941 mL
	10 mM	0.5947 mL	2.9735 mL	5.9471 mL
Please refer to the solubility information to select the appropriate solvent.				
In Vivo	<ol style="list-style-type: none"> <li>Add each solvent one by one: 10% DMSO &gt;&gt; 40% PEG300 &gt;&gt; 5% Tween-80 &gt;&gt; 45% saline Solubility: ≥ 2.5 mg/mL (14.87 mM); Clear solution</li> <li>Add each solvent one by one: 10% DMSO &gt;&gt; 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.5 mg/mL (14.87 mM); Clear solution</li> <li>Add each solvent one by one: 10% DMSO &gt;&gt; 90% corn oil Solubility: ≥ 2.5 mg/mL (14.87 mM); Clear solution</li> </ol>			

### BIOLOGICAL ACTIVITY

Description	5-Methoxysalicylic acid (5-MeOSA) is a natural compound, used as a useful matrix in the MALDI MS analysis of oligonucleotides when combined with spermine <sup>[1]</sup> .
IC <sub>50</sub> & Target	Human Endogenous Metabolite

### REFERENCES

---

[1]. Lee D, et al. 5-Methoxysalicylic acid matrix for ganglioside analysis with matrix-assisted laser desorption/ionization mass spectrometry. J Am Soc Mass Spectrom. 2015 Mar;26(3):522-5.

---

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: [tech@MedChemExpress.com](mailto:tech@MedChemExpress.com)

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA