



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## L-Thyroxine sodium xhydrate

<b>Cat. No.:</b>	HY-W010669		
<b>CAS No.:</b>	25416-65-3		
<b>Molecular Formula:</b>	C <sub>15</sub> H <sub>12</sub> I <sub>4</sub> NNaO <sub>5</sub>		
<b>Molecular Weight:</b>	816.87		
<b>Target:</b>	Thyroid Hormone Receptor		
<b>Pathway:</b>	Vitamin D Related/Nuclear Receptor		
<b>Storage:</b>	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	6 months
		-20°C	1 month

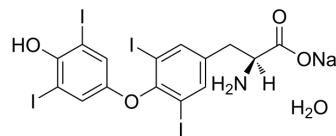

### BIOLOGICAL ACTIVITY

<b>Description</b>	L-Thyroxine sodium xhydrate (Levothyroxine; T4) is a synthetic hormone in the treatment of hypothyroidism. DIO enzymes convert biologically active hormone (Triiodothyronine, T3) from L-Thyroxine (T4).
<b>In Vivo</b>	<p>Deiodinases (DIOs), which catalyze the conversion of thyroxine (pro-hormone) to the active thyroid hormone, are associated with thyroid stimulating hormone (TSH) levels. DIO1 and DIO2 catalyze activation of thyroid hormone secretion in contrast to DIO3 playing role inactivation of the secretion. Activities of DIO1 and DIO2 play pivotal role in the negative feedback regulation of pituitary TSH secretion<sup>[1]</sup>. L-Thyroxine (T4) and Triiodothyronine (T3) hormones are known to modulate the expression of ionic channels, pumps and regulatory contractile proteins. Moreover, thyroid hormones have been shown to influence calcium homeostasis and flux responsible for excitation and contractility, with L-Thyroxine and Triiodothyronine modulating its pharmacological control and secretion. In rats fed 12 weeks with the iodine-free diet, a significant decrease in the levels of both Triiodothyronine and L-Thyroxine is observed when compared to the control group fed with standard diet (p&lt;0.001). In the group treated with low doses of L-Thyroxine, an increase in L-Thyroxine levels is observed (p=0.02) while Triiodothyronine levels remain virtually similar to the control group (p=0.19). Rats treated with high doses of L-Thyroxine display a significant increase in both Triiodothyronine and L-Thyroxine circulating concentrations compared to the non-treated hypothyroid group (p&lt;0.001 and p=0.004, respectively) and a significant increase in L-Thyroxine levels when compared to the control values (p=0.03)<sup>[2]</sup>.</p> <p>MCE has not independently confirmed the accuracy of these methods. They are for reference only.</p>

### PROTOCOL

<b>Animal Administration</b> <sup>[2]</sup>	<p>Rats<sup>[2]</sup></p> <p>Sprague-Dawley female rats (N=22) are used. Non-pregnant rats are divided into four groups: 1) control, 2) hypothyroidism, 3) hypothyroidism treated with low doses of L-Thyroxine (20 µg/kg/day) and 4) with high doses of L-Thyroxine (100 µg/kg/day). Control rats (group 1) are fed with standard diet, while the intervention rats are fed with iodine-free diet for 12 weeks to induce hypothyroidism (groups 2-4) which is continued for four more weeks to allow screening of hypothyroid status and L-Thyroxine-treatment. Food and water (iodine-free diet) are available ad libitum. The hypothyroid group treated with low (group 3) or high doses of L-Thyroxine (group 4) are injected intraperitoneally every 24 h with respectively 20 µg/kg/day and 100 µg/kg/day. Blood samples are collected for thyroid function screening at week 12 and 16 following the initiation of either the control or iodine-free diet. Hysterectomy is performed under general anesthesia (isoflurane 2%) at</p>
---	---

---

the end of the treatment and the two uterine horns are placed in physiological Krebs' solution until isometric tension measurements within no more than 1 h.  
MCE has not independently confirmed the accuracy of these methods. They are for reference only.

---

## REFERENCES

---

- [1]. Arici M, et al. Association between genetic polymorphism and levothyroxine bioavailability in hypothyroid patients. *Endocr J*. 2018 Jan 11.
- [2]. Corriveau S, et al. Levothyroxine treatment generates an abnormal uterine contractility patterns in an in vitro animal model. *J Clin Transl Endocrinol*. 2015 Sep 9;2(4):144-149.
- 

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: [tech@MedChemExpress.com](mailto:tech@MedChemExpress.com)

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA