



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

2-Hydroxy-6-methoxybenzoic acid

Cat. No.:	HY-W017100		
CAS No.:	3147-64-6		
Molecular Formula:	C ₈ H ₈ O ₄		
Molecular Weight:	168.15		
Target:	Others		
Pathway:	Others		
Storage:	Powder	-20°C	3 years
		4°C	2 years
	In solvent	-80°C	2 years
		-20°C	1 year

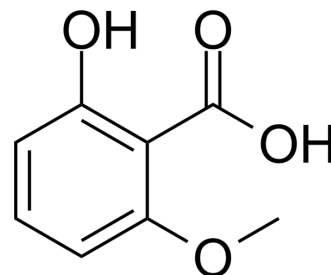

SOLVENT & SOLUBILITY

In Vitro	DMSO : 100 mg/mL (594.71 mM; Need ultrasonic)			
		Solvent Concentration	Mass	
			1 mg	5 mg
			10 mg	
Preparing Stock Solutions	1 mM	5.9471 mL	29.7354 mL	59.4707 mL
	5 mM	1.1894 mL	5.9471 mL	11.8941 mL
	10 mM	0.5947 mL	2.9735 mL	5.9471 mL
Please refer to the solubility information to select the appropriate solvent.				
In Vivo	<ol style="list-style-type: none"> Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline Solubility: ≥ 2.5 mg/mL (14.87 mM); Clear solution Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline) Solubility: ≥ 2.5 mg/mL (14.87 mM); Clear solution Add each solvent one by one: 10% DMSO >> 90% corn oil Solubility: ≥ 2.5 mg/mL (14.87 mM); Clear solution 			

BIOLOGICAL ACTIVITY

Description	2-Hydroxy-6-methoxybenzoic acid can be used for the determination of acetylsalicylic acid and its major metabolite, salicylic acid, in animal plasma. 2-Hydroxy-6-methoxybenzoic acid exhibits significant analgesic effects ^{[1][2]} .
--------------------	--

REFERENCES

[1]. Xu X, et al. Rapid and sensitive determination of acetylsalicylic acid and salicylic acid in plasma using liquid chromatography-tandem mass spectrometry: application

to pharmacokinetic study. Biomed Chromatogr. 2009 Sep;23(9):973-9.

[2]. Palacios-Espinosa F, et al. Antinociceptive, hypoglycemic and spasmolytic effects of Brickellia veronicifolia. J Ethnopharmacol. 2008 Aug 13;118(3):448-54.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: tech@MedChemExpress.com

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA