



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

## Ethylene glycol-d<sub>4</sub>

<b>Cat. No.:</b>	HY-Y0338S2	
<b>CAS No.:</b>	2219-51-4	
<b>Molecular Formula:</b>	C <sub>2</sub> H <sub>2</sub> D <sub>4</sub> O <sub>2</sub>	
<b>Molecular Weight:</b>	66.09	
<b>Target:</b>	Endogenous Metabolite	
<b>Pathway:</b>	Metabolic Enzyme/Protease	
<b>Storage:</b>	Pure form	-20°C 3 years 4°C 2 years
	In solvent	-80°C 6 months -20°C 1 month

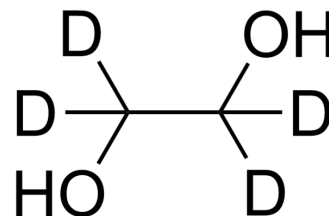

### SOLVENT & SOLUBILITY

#### In Vitro

H<sub>2</sub>O : 200 mg/mL (3026.18 mM; Need ultrasonic)  
DMSO : 100 mg/mL (1513.09 mM; Need ultrasonic)

	Solvent		Mass		
	Concentration		1 mg	5 mg	10 mg
Preparing Stock Solutions	1 mM		15.1309 mL	75.6544 mL	151.3088 mL
	5 mM		3.0262 mL	15.1309 mL	30.2618 mL
	10 mM		1.5131 mL	7.5654 mL	15.1309 mL

Please refer to the solubility information to select the appropriate solvent.

#### In Vivo

- Add each solvent one by one: 10% DMSO >> 40% PEG300 >> 5% Tween-80 >> 45% saline  
Solubility: ≥ 2.5 mg/mL (37.83 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% (20% SBE-β-CD in saline)  
Solubility: ≥ 2.5 mg/mL (37.83 mM); Clear solution
- Add each solvent one by one: 10% DMSO >> 90% corn oil  
Solubility: ≥ 2.5 mg/mL (37.83 mM); Clear solution

### BIOLOGICAL ACTIVITY

#### Description

Ethylene glycol-d<sub>4</sub> is the deuterium labeled Ethylene glycol[1].

#### In Vitro

Stable heavy isotopes of hydrogen, carbon, and other elements have been incorporated into drug molecules, largely as tracers for quantitation during the drug development process. Deuteration has gained attention because of its potential to affect the pharmacokinetic and metabolic profiles of drugs<sup>[1]</sup>.

MCE has not independently confirmed the accuracy of these methods. They are for reference only.

---

## REFERENCES

---

[1]. Russak EM, et al. Impact of Deuterium Substitution on the Pharmacokinetics of Pharmaceuticals. *Ann Pharmacother.* 2019 Feb;53(2):211-216.

---

**Caution: Product has not been fully validated for medical applications. For research use only.**

Tel: 609-228-6898

Fax: 609-228-5909

E-mail: [tech@MedChemExpress.com](mailto:tech@MedChemExpress.com)

Address: 1 Deer Park Dr, Suite Q, Monmouth Junction, NJ 08852, USA