

# Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



## Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

## Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

## SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in





www.rockland.com tech@rockland.com +1 484.791.3823

#### Datasheet for 009-0033

## **Human Albumin Rhodamine**

#### **Overview**

| Description:  | Human Albumin Rhodamine Conjugated - 009-0033 |
|---------------|---|
| Item No.:     | 009-0033                                      |
| Size:         | 1 mg  |
| Applications: | Dot Blot                                      |
| Origin:       | Human   |

### **Product Details**

| Background:        | Human albumin or serum albumin is encoded by the ALB gene and is the most abundant plasma protein in mammals. Human albumin is essential for maintaining the osmotic pressure needed for proper distribution of body fluids between intravascular compartments and body tissues. Human albumin also acts as a plasma carrier by non-specifically binding several hydrophobic steroid hormones and as a transport protein for hemin and fatty acids. Too much serum albumin in the body can be harmful. |
|--------------------|--|
| Synonyms:          | Human Albumin rhodamine conjugation, Human Albumin TRITC conjugation   |
| Species of Origin: | Human  |
| Conjugate:         | Rhodamine (TRITC)  |
| Format:            | Albumin  |
| Type:              | Native Protein   |
| F/P Ratio:         | 1.6  |
| Specific Activity: | 1.2  |
|                    |  |

### **Target Details**

Purity/Specificity: Human Albumin Rhodamine was prepared from chromatographically pure Human Albumin.

Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Human Albumin.

## **Application Details**

www.rockland.com Page 1 of 3





www.rockland.com tech@rockland.com +1 484.791.3823

| <b>Tested Applications:</b> | Dot Blot  |
|-----------------------------|---|
| Application Note:           | Human Albumin Rhodamine has been tested by dot blot and is designed for immunofluorescence microscopy, fluorescence based plate assays (FLISA), fluorescent western blotting, multiplex analysis, including multicolor imaging, utilizing various commercial platforms. |
| Assay Dilutions:            | All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.   |

### **Formulation**

| Physical State:        | Lyophilized  |
|------------------------|--|
| Concentration:         | 1.0 mg/mL by UV absorbance at 280 nm                       |
| Buffer:                | 0.02 M Potassium Phosphate, 0.15 M Sodium Chloride, pH 7.2 |
| Preservative:          | 0.01% (w/v) Sodium Azide                                   |
| Stabilizer:            | 10 mg/ml Polyethylene Glycol (PEG-8000)                    |
| Reconstitution Volume: | 1.0 mL   |
| Reconstitution Buffer: | Restore with deionized water (or equivalent)               |

## **Shipping & Handling**


| <b>Shipping Condition:</b> | Ambient   |
|----------------------------|---|
| Storage Condition:         | Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use. |
| Expiration:                | Expiration date is one (1) year from date of receipt.   |

## **Images**

www.rockland.com Page 2 of 3




www.rockland.com tech@rockland.com +1 484.791.3823

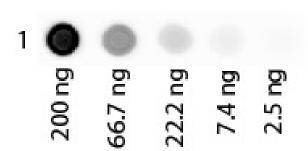

#### **Dot Blot**

Dot Blot of Human Albumin Rhodamine Conjugated. Antigen: Rhodamine Conjugated Human Albumin. Load: 3-fold serial dilution starting at 200 ng. Primary antibody: None. Secondary antibody: None. Block: MB-070 for 1 H at RT.

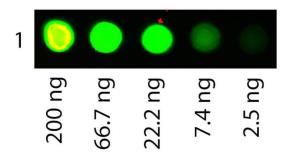

#### **Dot Blot**

Dot Blot of Human Albumin Rhodamine Conjugated. Antigen: Rhodamine Conjugated Human Albumin. Load: 3-fold serial dilution starting at 200 ng. Primary antibody: None. Secondary antibody: None. Block: MB-070 for 1 H at RT.

#### **Disclaimer**

No test method can provide total assurance that the hepatitis B virus, hepatitis C virus, human immunodeficiency virus, or any other infectious agents are absent. Thus, all blood products, including purified proteins derived from human blood sources, should be handled at Biosafety Level 2 as recommended by the CDC\NIH manual entitled Biosafety in Microbiological and Biomedical Laboratories for potentially infectious human serum, blood specimens or proteins derived from same. Source material for the human blood product supplied to your facility has been tested for the detection of HIV antibody, Hepatitis B surface antigen, antibody to Hepatitis C, HIV 1 antigen(s), antibody to HTLV - I/II, and syphilis by FDA guidelines. All units were found to be non-reactive/negative for these tests. All human blood source material is collected in FDA licensed centers and is tested with FDA approved test kits.

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.

www.rockland.com Page 3 of 3