



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Datasheet for 009-0134

Human Transferrin

Overview

Description:	Human Transferrin- 009-0134
Item No.:	009-0134
Size:	10 mg
Applications:	SDS-PAGE
Origin:	Human

Product Details

Background:	Human transferrin is encoded by the TF gene and is an iron-binding blood plasma glycoprotein that controls the level of free iron in biological fluids. Human transferrin binds iron very tightly but reversibly. Human transferrin is the most important iron pool in mammals. Human transferrin has a molecular weight of around 80 kDa and contains 2 specific high-affinity Fe(III) binding sites. The affinity of Human transferrin for Fe(III) is extremely high but decreases progressively with decreasing pH below neutrality.
Synonyms:	apotransferrin, Transferrin Control, Siderophilin, TRF, apo-hTF
Species of Origin:	Human
Format:	Transferrin
Type:	Native Protein

Target Details

Purity/Specificity:	Human transferrin was prepared from normal serum by a multi-step process which includes delipidation and selective precipitation followed by extensive dialysis. Assay by immunoelectrophoresis resulted in a single precipitin arc against anti-Human Transferrin and anti-Human Serum.
Relevant Links:	<ul style="list-style-type: none">• NCBI - AAB22049.1• UniProtKB - Q06AH7• GenelD - 7018

Application Details

Tested Applications:	SDS-PAGE
Application Note:	Human Apo Transferrin from human serum was lyophilized from 0.2M ammonium bicarbonate for use in cell culture and assays. Transferrin has been tested in SDS-Page and is also suitable for use as antigen or ligand in immunochemical reactions, as a control or standard in assays, for conjugation and most other immunological methods requiring highly purified proteins.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	1:10,000 - 1:20,000
WB:	1:2,000-1:10,000

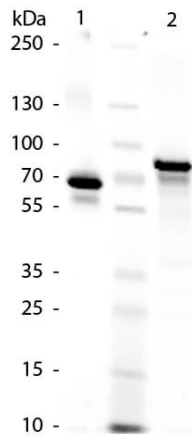
Formulation

Physical State:	Lyophilized
Concentration:	10.0 mg/mL by UV absorbance at 280 nm
Buffer:	See application note.
Preservative:	None
Stabilizer:	None
Reconstitution Volume:	1.0 mL
Reconstitution Buffer:	Restore with deionized water (or equivalent)

Shipping & Handling

Shipping Condition:	Ambient
Storage Condition:	Store vial at 4° C prior to restoration. For extended storage aliquot contents and freeze at -20° C or below. Avoid cycles of freezing and thawing. Centrifuge product if not completely clear after standing at room temperature. This product is stable for several weeks at 4° C as an undiluted liquid. Dilute only prior to immediate use.
Expiration:	Expiration date is one (1) year from date of receipt.

Images

**SDS-PAGE**

SDS-Page of Human Transferrin. Lane 1: Human Transferrin – Reduced. Lane 2: Human Transferrin – Non-Reduced.
Load: 1.0 µg per lane. Observed/Predicted Size: ~70 kDa for Reduced, ~80 kDa for Non-Reduced Transferrin.

Disclaimer

No test method can provide total assurance that the hepatitis B virus, hepatitis C virus, human immunodeficiency virus, or any other infectious agents are absent. Thus, all blood products, including purified proteins derived from human blood sources, should be handled at Biosafety Level 2 as recommended by the CDC\NIH manual entitled Biosafety in Microbiological and Biomedical Laboratories for potentially infectious human serum, blood specimens or proteins derived from same. Source material for the human blood product supplied to your facility has been tested for the detection of HIV antibody, Hepatitis B surface antigen, antibody to Hepatitis C, HIV 1 antigen(s), antibody to HTLV - I/II, and syphilis by FDA guidelines. All units were found to be non-reactive/negative for these tests. All human blood source material is collected in FDA licensed centers and is tested with FDA approved test kits.

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.