



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Datasheet for TMBE-100

TMBE ELISA Peroxidase Substrate**Overview**

Description:	TMBE ELISA Peroxidase Substrate - TMBE-100
Item No.:	TMBE-100
Size:	100 mL
Applications:	ELISA

Product Details

Background:	TMBE (3,3',5,5'-Tetramethylbenzidine) is a colorimetric substrate for peroxidase (HRP)-based enzyme immunoassays, notably ELISAs, and is known for producing a blue color that changes to yellow upon stopping the reaction with sulfuric acid. The color change is quantifiable at 450 nm. TMBE's high sensitivity makes it ideal for detecting low concentrations of analytes.
Synonyms:	3,3',5,5'-Tetramethylbenzidine, TMBE ELISA substrate, TMBE Substrate, chromogenic ELISA substrate

Target Details

Purity/Specificity:	pH: 4.0 +/- 0.2 Stability @ 18° to 26°C: PASS Stability @ 4°C: PASS QC Raw Material: PASS Absorbance check of final product: PASS Performance Data per ELISA: PASS
Relevant Links:	<ul style="list-style-type: none">TMBE SDS

Application Details

Tested Applications:	ELISA
-----------------------------	-------

Application Note:	TMB ELISA Peroxidase Substrate comes ready to use. No dilutions are required. TMB ELISA Peroxidase Substrate will produce a soluble blue end product read at 370 nm or 655 nm. TMB ELISA Peroxidase Substrate incubation time will vary depending on the assay conditions.
Assay Dilutions:	All assays should be optimized by the user. Recommended dilutions (if any) may be listed below.
ELISA:	1X


Formulation

Physical State:	Liquid - clear to very light blue colored liquid
Concentration:	1x
Stabilizer:	Proprietary buffer with enhancer

Shipping & Handling

Shipping Condition:	Ambient
Storage Condition:	Store container at 4° C prior to opening. Protect from moisture and light. No special shipping conditions or precautions are required.
Expiration:	Expiration date is one (1) year from date of receipt.

Images

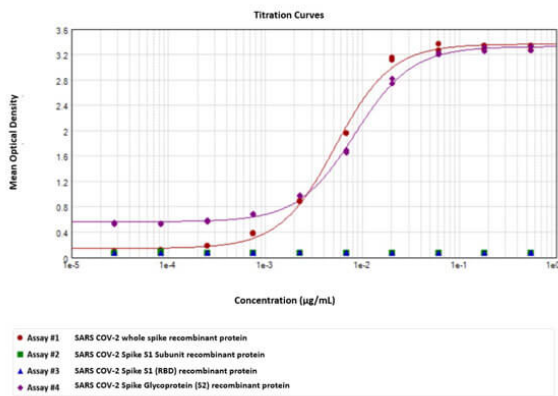


ELISA

ELISA results of Mouse Anti-SARS CoV-2 Spike Protein (S2) Antibody using TMBE-1000. Each well was coated in duplicate with 0.5µg of rec SARS CoV-2 Whole Spike protein [RED line], rec SARS CoV-2 Spike (S1) [GREEN line], rec SARS CoV-2 Spike (S1) RBD protein [BLUE line], and rec SARS CoV-2 Spike (S2) glycoprotein [PURPLE line]. The working dilution for SARS-CoV-2 whole spike is 1:190,000 and SARS CoV-2 Spike (S2) glycoprotein is 1:120,000. The starting dilution of antibody was 0.55µg/ml and the X-axis represents the Log10 of a 3-fold dilution. This titration is a 4-parameter curve fit where the IC50 is defined as the titer of the antibody. Assay performed using Rabbit Anti-Mouse IgG HRP conjugated (p/n 610-403-C46) at 1:8000 and TMB substrate (p/n TMBE-1000). SARS CoV-2 Spike (S2) EC50: 8ng/mL; SARS CoV-2 Whole Spike EC50: 5ng/mL.

**Bottle**

TMBS ELISA Peroxidase Substrate

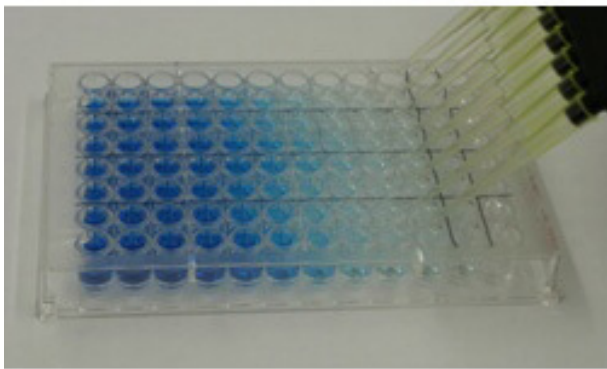**ELISA**

Rockland Immunochemicals produces a wide variety of buffers and substrates for use in ELISAs. Antigen was diluted in ELISA Microwell Coating Stabilizer (p/n MB-063-0100) added to the microwell plate and incubated overnight at 4°C. The plate was then blocked with ELISA Microwell Blocking Buffer with Stabilizer (p/n MB-064-1000) for 2 hours. The primary antibody was diluted in PBS Fish Gel Concentrate (1:10)(p/n MB-066-0100), added to the plate, and allowed to incubate 1 hour at room temperature. HRP conjugated secondary antibody was diluted in HRP Conjugate Stabilizer (p/n MB-060-0100), added to the plate, and allowed to incubate for 30 minutes at room temperature. TMB ELISA Peroxidase Substrate (p/n TMBS-1000) was added to the plate and allowed to incubate for 30 minutes at room temperature. The reaction was then stopped with 1M HCl and read at 450nm.

References

- Escalera A et al. SARS-CoV-2 infection induces robust mucosal antibody responses in the upper respiratory tract. *iScience*. (2024)
- Kukla, DA et al. Assessing the compatibility of primary human hepatocyte culture within porous silk sponges. *Rsc Advances* (2022)
- Quach HQ et al. Detection of SARS-CoV-2 peptide-specific antibodies in Syrian hamster serum by ELISA. *J Immunol Methods*. (2022)
- Pandit M et al. AMPK suppresses Th2 cell responses by repressing mTORC2. *Exp Mol Med*. (2022)
- Pandit M et al. LKB1-PTEN axis controls Th1 and Th17 cell differentiation via regulating mTORC1. *J Mol Med (Berl)*. (2021)

- Kratzer B et al. Immunological imprint of COVID-19 on human peripheral blood leukocyte populations. *Allergy* (2021)
- Scharer CD et al. Antibody-secreting cell destiny emerges during the initial stages of B-cell activation. *Nat Commun.* (2020)
- Shao W et al. Serum lipoprotein–derived fatty acids regulate hypoxia-inducible factor. *J Biol Chem.* (2020)
- Harshbarger W et al. Convergent structural features of respiratory syncytial virus neutralizing antibodies and plasticity of the site V epitope on prefusion F. *PLoS Pathog.* (2020)
- Kukla DA. et. al. Microscale Collagen and Fibroblast Interactions Enhance Primary Human Hepatocyte Functions in Three-Dimensional Models. *Gene Expr.* (2020)
- Kukla DA. et. al. Primary Human Hepatocytes Maintain Long-term Functions in Porous Silk Scaffolds Containing Extracellular Matrix Proteins. *bioRxiv* (2020)
- Ferdman J et al. Intra-seasonal antibody repertoire analysis of a subject immunized with an MF59®-adjuvanted pandemic 2009 H1N1 vaccine. *Vaccine.* (2019)
- Haines, RR et al. The Histone Demethylase LSD1 Regulates B Cell Proliferation and Plasmablast Differentiation. *Journal of Immunology (Baltimore, Md. : 1950)* (2018)
- Tomita, K et al. Antibodies against adenovirus fiber and penton base proteins inhibit adenovirus vector-mediated transduction in the liver following systemic administration. *Scientific Reports* (2018)
- Barwick, BG et al. B cell activation and plasma cell differentiation are inhibited by de novo DNA methylation. *Nature Communications* (2018)
- Steff, AM et al. Pre-fusion RSV F strongly boosts pre-fusion specific neutralizing responses in cattle pre-exposed to bovine RSV. *Nature Communications* (2017)
- O'shannessy, DJ et al. Novel antibody probes for the characterization of endosialin/TEM-1. *Oncotarget* (2016)
- Parry, HM et al. Cytomegalovirus viral load within blood increases markedly in healthy people over the age of 70 years. *Immunity & Ageing : I & A* (2016)
- Ciferri et al. Antigenic Characterization of the HCMV gH/gL/gO and Pentamer Cell Entry Complexes Reveals Binding Sites for Potently Neutralizing Human Antibodies. *PLOS Pathogens* (2015)
- Shimizu, K et al. Suppression of leaky expression of adenovirus genes by insertion of microRNA-targeted sequences in the replication-incompetent adenovirus vector genome. *Molecular Therapy. Methods & Clinical Development* (2014)
- O'shannessy, DJ et al. Characterization of the human folate receptor alpha via novel antibody-based probes. *Oncotarget* (2011)

Disclaimer

This product is for research use only and is not intended for therapeutic or diagnostic applications. Please contact a technical service representative for more information. All products of animal origin manufactured by Rockland Immunochemicals are derived from starting materials of North American origin. Collection was performed in United States Department of Agriculture (USDA) inspected facilities and all materials have been inspected and certified to be free of disease and suitable for exportation. All properties listed are typical characteristics and are not specifications. All suggestions and data are offered in good faith but without guarantee as conditions and methods of use of our products are beyond our control. All claims must be made within 30 days following the date of delivery. The prospective user must determine the suitability of our materials before adopting them on a commercial scale. Suggested uses of our products are not recommendations to use our products in violation of any patent or as a license under any patent of Rockland Immunochemicals, Inc. If you require a commercial license to use this material and do not have one, then return this material, unopened to: Rockland Inc., P.O. BOX 5199, Limerick, Pennsylvania, USA.