



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

Calpastatin (m): 293T Lysate: sc-118967

BACKGROUND

Calpains are nonlysosomal, calcium-activated intracellular cysteine proteases that mediate specific Ca^{2+} -dependent processes including cell fusion, mitosis and meiosis. Calpains are heterodimers of a small regulatory subunit and one of three large catalytic subunits, designated Calpain 1, Calpain 2 and Calpain p94. Calpain 1 is an intracellular calcium-dependent protease that cleaves cytoskeletal and submembranous proteins. Calpain 1 co-localizes with human leukocyte antigen-DR (HLA-DR) on activated microglia in the aging brain. Calpain influences the process of spermatogenesis and the events preceding fertilization, such as the acrosome reaction. Calpastatin regulates Calpain by inhibiting both the proteolytic activity of Calpain and its binding to membranes. Calpastatin exists in two types, tissue type and erythrocyte type, resulting from both alternative splicing and proteolytic processing.

REFERENCES

- Murachi, T. 1984. Calcium-dependent proteinases and specific inhibitors: Calpain and Calpastatin. *Biochem. Soc. Symp.* 45: 149-167.
- Takano, E., Ueda, M., Tsunekawa, S., Murakami, T., Maki, M., Hatanaka, M. and Murachi, T. 1991. Molecular diversity of erythrocyte Calpastatin. *Biomed. Biochim. Acta* 50: 517-521.
- Takano, E., Nosaka, T., Lee, W.J., Nakamura, K., Takahashi, T., Funaki, M., Okada, H., Hatanaka, M. and Maki, M. 1993. Molecular diversity of Calpastatin in human erythroid cells. *Arch. Biochem. Biophys.* 303: 349-354.
- Kawasaki, H. and Kawashima, S. 1996. Regulation of the Calpain-Calpastatin system by membranes. *Mol. Mem. Biol.* 13: 217-224.
- Johnson, G.V. and Guttman, R.P. 1997. Calpains: intact and active? *Bioessays* 19: 1011-1018.
- Elce, J.S., Hegadorn, C. and Arther, J.S.C. 1997. Autolysis, Ca^{2+} requirement, and heterodimer stability in μ -Calpain. *J. Biol. Chem.* 272: 11268-11275.
- Niapour, M., Yu, Y. and Berger, S.A. 2008. Regulation of Calpain activity by c-Myc through Calpastatin and promotion of transformation in c-Myc-negative cells by Calpastatin suppression. *J. Biol. Chem.* 283: 21371-21381.
- Kiss, R., Kovács, D., Tompa, P. and Perczel, A. 2008. Local structural preferences of Calpastatin, the intrinsically unstructured protein inhibitor of Calpain. *Biochemistry.* 47: 6936-45.
- Kiss, R., Bozoky, Z., Kovács, D., Róna, G., Friedrich, P., Dvortsák, P., Weisemann, R., Tompa, P. and Perczel, A. 2008. Calcium-induced tripartite binding of intrinsically disordered calpastatin to its cognate enzyme, calpain. *FEBS Lett.* 582: 2149-54.

CHROMOSOMAL LOCATION

Genetic locus: Cast (mouse) mapping to 13 C1.

PRODUCT

Calpastatin (m): 293T Lysate represents a lysate of mouse Calpastatin transfected 293T cells and is provided as 100 μ g protein in 200 μ l SDS-PAGE buffer.

APPLICATIONS

Calpastatin (m): 293T Lysate is suitable as a Western Blotting positive control for mouse reactive Calpastatin antibodies. Recommended use: 10-20 μ l per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

Calpastatin (PI-11): sc-32324 is recommended as a positive control antibody for Western Blot analysis of enhanced mouse Calpastatin expression in Calpastatin transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).


RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:
1) Western Blotting: use m-IgG κ BP-HRP: sc-516102 or m-IgG κ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

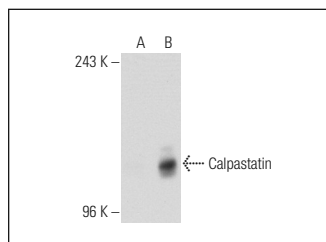

Calpastatin (PI-11): sc-32324. Western blot analysis of Calpastatin expression in non-transfected: sc-117752 (A) and mouse Calpastatin transfected: sc-118967 (B) 293T whole cell lysates.

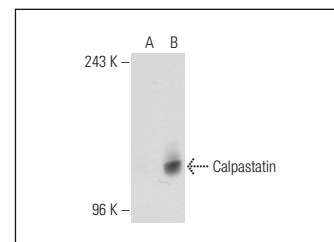

Calpastatin (2Q31): sc-70486. Western blot analysis of Calpastatin expression in non-transfected: sc-117752 (A) and mouse Calpastatin transfected: sc-118967 (B) 293T whole cell lysates.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.