



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

HSF1 (m): 293T Lysate: sc-120904

BACKGROUND

Prokaryotic and eukaryotic cells respond to thermal and chemical stress by inducing a group of genes collectively designated heat shock genes. In eukaryotes, this gene expression is regulated primarily at the transcription level. Heat shock transcription factors 1 and 2, also designated HSF1 and 2, are involved in this regulation. HSF1 and HSF2 are upregulated by estrogen at both the mRNA and protein level. HSF1 is normally found as a monomer, whose transcriptional activity is repressed by constitutive phosphorylation. Upon activation, HSF1 forms trimers, gains DNA binding activity and is translocated to the nucleus. HSF2 activity is associated with differentiation and development, and, like HSF1, binds DNA as a trimer. Both HSF1 and HSF2 are known to be induced by proteasome inhibitors of the ubiquitin pathway.

REFERENCES

1. Tanguay, R.M. 1988. Transcriptional activation of heat shock genes in eukaryotes. *Biochem. Cell Biol.* 66: 584-593.
2. Yang, X., et al. 1995. Estrogen dependent expression of heat shock transcription factor: implications for uterine synthesis of heat shock proteins. *J. Steroid Biochem. Mol. Biol.* 52: 415-419.
3. Wyman, C., et al. 1995. Determination of HSF2 stoichiometry at looped DNA complexes using scanning force microscopy. *EMBO J.* 14: 117-123.
4. Rallu, M., et al. 1997. Function and regulation of HSF2 during mouse embryogenesis. *Proc. Natl. Acad. Sci. USA* 94: 2392-2397.
5. He, B., et al. 1998. Glycogen synthase kinase 3b and extracellular signal-regulated kinase inactivate HSF1 by facilitating the disappearance of transcriptionally active granules after heat shock. *Mol. Cell Biol.* 18: 6624-6633.
6. Kawazoe, Y., et al. 1998. Proteasome inhibition leads to the activation of all members of the heat shock factor family. *Eur. J. Biochem.* 255: 356-362.
7. Mathew, A., et al. 1998. Heat shock response and protein degradation: regulation of HSF2 by the ubiquitin-proteasome pathway. *Mol. Cell Biol.* 18: 5091-5098.

CHROMOSOMAL LOCATION

Genetic locus: Hsf1 (mouse) mapping to 15 D3.

PRODUCT

HSF1 (m): 293T Lysate represents a lysate of mouse HSF1 transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

APPLICATIONS

HSF1 (m): 293T Lysate is suitable as a Western Blotting positive control for mouse reactive HSF1 antibodies. Recommended use: 10-20 µl per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

HSF1 (C-5): sc-17756 is recommended as a positive control antibody for Western Blot analysis of enhanced mouse HSF1 expression in HSF1 transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:
 1) Western Blotting: use m-IgGκ BP-HRP: sc-516102 or m-IgGκ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

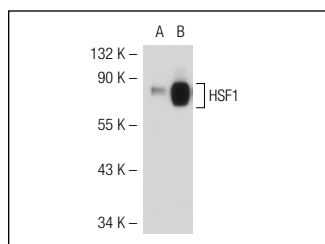

HSF1 (C-5): sc-17756. Western blot analysis of HSF1 expression in non-transfected: sc-117752 (A) and mouse HSF1 transfected: sc-120904 (B) 293T whole cell lysates.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.