



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

MMP-12 (m): 293T Lysate: sc-121694

BACKGROUND

The matrix metalloproteinases (MMP) are a family of peptidase enzymes responsible for the degradation of extracellular matrix components, including collagen, gelatin, Fibronectin, Laminin and proteoglycan. Transcription of MMP genes is differentially activated by phorbol ester, lipopolysaccharide (LPS) or staphylococcal enterotoxin B (SEB). MMP catalysis requires both calcium and zinc. MMP-12 (also designated macrophage metalloelastase) is produced in alveolar macrophages and degrades elastin. MMP-12 may contribute to elastin degradation occurring in granulomatous skin diseases and may also participate in macrophage migration through the epidermal and vascular basement membranes in inflammatory disorders.

REFERENCES

1. Shapiro, S.D., Griffin, G.L., Gilbert, D.J., Jenkins, N.A., Copeland, N.G., Welgus, H.G., Senior, R.M. and Ley, T.J. 1992. Molecular cloning, chromosomal localization, and bacterial expression of a murine macrophage metalloelastase. *J. Biol. Chem.* 267: 4664-4671.
2. Birkedal-Hansen, H., Moore, W.G., Bodden, M.K., Windsor, L.J., Birkedal-Hansen, B., DeCarlo, A. and Engler, J.A. 1993. Matrix metalloproteinases: a review. *Crit. Rev. Oral Biol. Med.* 4: 197-250.
3. Shapiro, S.D., Kobayashi, D.K. and Ley, T.J. 1993. Cloning and characterization of a unique elastolytic metalloproteinase produced by human alveolar macrophages. *J. Biol. Chem.* 268: 23824-23829.
4. Reinemer, P., Grams, F., Huber, R., Kleine, T., Schnierer, S., Piper, M., Tschesche, H. and Bode, W. 1994. Structural implications for the role of the N-terminus in the "superactivation" of collagenases. A crystallographic study. *FEBS Lett.* 338: 227-233.
5. Machein, U. and Conca, W. 1997. Expression of several matrix metalloproteinase genes in human monocytic cells. *Adv. Exp. Med. Biol.* 421: 247-251.
6. Vaalamo, M., Kariniemi, A.L., Shapiro, S.D. and Saarialho-Kere, U. 1999. Enhanced expression of human metalloelastase (MMP-12) in cutaneous granulomas and macrophage migration. *J. Invest. Dermatol.* 112: 499-505.

CHROMOSOMAL LOCATION

Genetic locus: *Mmp12* (mouse) mapping to 9 A1.

PRODUCT

MMP-12 (m): 293T Lysate represents a lysate of mouse MMP-12 transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

APPLICATIONS

MMP-12 (m): 293T Lysate is suitable as a Western Blotting positive control for mouse reactive MMP-12 antibodies. Recommended use: 10-20 µl per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

MMP-12 (H-9): sc-390284 is recommended as a positive control antibody for Western Blot analysis of enhanced mouse MMP-12 expression in MMP-12 transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:
 1) Western Blotting: use m-IgGκ BP-HRP: sc-516102 or m-IgGκ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

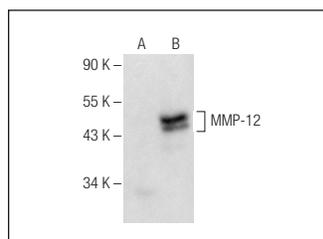

MMP-12 (H-9): sc-390284. Western blot analysis of MMP-12 expression in non-transfected: sc-117752 (A) and mouse MMP-12 transfected: sc-121694 (B) 293T whole cell lysates.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.