



**SZABO  
SCANDIC**

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](http://linkedin.com/company/szaboscandic)



# SSSCA1 (m): 293T Lysate: sc-123793

## BACKGROUND

SSSCA1 (Sjogren syndrome/scleroderma autoantigen 1), also known as p27, is a 199 amino acid protein that may play a role in mitosis. Considered a centromere-associated protein, SSSCA1 may induce anti-centromere antibodies. SSSCA1 is phosphorylated upon DNA damage by ATM or ATR. The gene encoding SSSCA1 maps to human chromosome 11, which houses over 1,400 genes and comprises nearly 4% of the human genome. Jervell and Lange-Nielsen syndrome, Jacobsen syndrome, Niemann-Pick disease, hereditary angioedema and Smith-Lemli-Opitz syndrome are associated with defects in genes that map to chromosome 11.

## REFERENCES

1. Grossfeld, P.D., et al. 2004. The 11q terminal deletion disorder: a prospective study of 110 cases. *Am. J. Med. Genet. A* 129: 51-61.
2. Loussouarn, G., et al. 2006. KCNQ1 K<sup>+</sup> channel-mediated cardiac channelopathies. *Methods Mol. Biol.* 337: 167-183.
3. Taylor, T.D., et al. 2006. Human chromosome 11 DNA sequence and analysis including novel gene identification. *Nature* 440: 497-500.
4. Zehlein, J., et al. 2006. Skipping of Exon 1 in the KCNQ1 gene causes Jervell and Lange-Nielsen syndrome. *J. Biol. Chem.* 281: 35397-35403.
5. Ataga, K.I., et al. 2007. β-thalassaemia and sickle cell anaemia as paradigms of hypercoagulability. *Br. J. Haematol.* 139: 3-13.
6. Berger, A.C., et al. 2007. The subcellular localization of the Niemann-Pick Type C proteins depends on the adaptor complex AP-3. *J. Cell Sci.* 120: 3640-3652.
7. Lee, J.H. and Paull, T.T. 2007. Activation and regulation of ATM kinase activity in response to DNA double-strand breaks. *Oncogene* 26: 7741-7748.
8. O'Connor, M.J., et al. 2007. Targeted cancer therapies based on the inhibition of DNA strand break repair. *Oncogene* 26: 7816-7824.
9. Kaste, S.C., et al. 2008. Wilms tumour: prognostic factors, staging, therapy and late effects. *Pediatr. Radiol.* 38: 2-17.

## CHROMOSOMAL LOCATION

Genetic locus: Sssc1 (mouse) mapping to 19 A.

## PRODUCT

SSSCA1 (m): 293T Lysate represents a lysate of mouse SSSCA1 transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

## APPLICATIONS

SSSCA1 (m): 293T Lysate is suitable as a Western Blotting positive control for mouse reactive SSSCA1 antibodies.

Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.


SSSCA1 (G-12): sc-515430 is recommended as a positive control antibody for Western Blot analysis of enhanced mouse SSSCA1 expression in SSSCA1 transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

## RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:

1) Western Blotting: use m-IgG<sub>X</sub> BP-HRP: sc-516102 or m-IgG<sub>X</sub> BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

## DATA

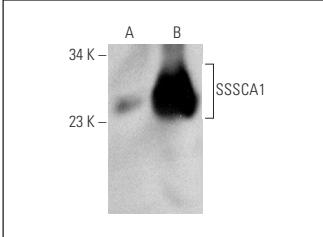

SSSCA1 (G-12): sc-515430. Western blot analysis of SSSCA1 expression in non-transfected: sc-117752 (**A**) and mouse SSSCA1 transfected: sc-123793 (**B**) 293T whole cell lysates.

## STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

## RESEARCH USE

For research use only, not for use in diagnostic procedures.

## PROTOCOLS

See our web site at [www.scbt.com](http://www.scbt.com) for detailed protocols and support products.