



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

EPI64 (m2): 293T Lysate: sc-125304

BACKGROUND

The Na⁺/H⁺ exchange protein (NHE) functions in transepithelial Na⁺ absorption and is primarily expressed in the intestinal and renal brush border membrane. NHE regulatory factor 1 (NHERF-1) interacts with NHE through two PDZ (for PSD-95, discs-large, and ZO-1 homology) domains, which are protein-protein interaction modules that associate with specific C-terminal motifs on target proteins. Also known as EBP50, NHERF-1 facilitates cAMP inhibition of NHE to decrease Na⁺ adsorption. NHERF-2, also known as E3KARP, is ubiquitously expressed as a protein which also functions in NHE-2 regulation. EBP-PDZ interactor (EPI64) contains a C-terminal -DTYL sequence that binds to the first PDZ domain of NHERF-1 and NHERF-2. EPI64 is ubiquitously expressed and localizes with NHERF-1 *in vitro*. The gene encoding human EPI64 maps to chromosome 22q12.2.

REFERENCES

- Sheng, M. 1996. PDZs and receptor/channel clustering: rounding up the latest suspects. *Neuron* 17: 575-578.
- Yun, C.H., Oh, S., Zizak, M., Steplock, D., Tsao, S., Tse, C.M., Weinman, E.J. and Donowitz, M. 1997. cAMP-mediated inhibition of the epithelial brush border Na⁺/H⁺ exchanger, NHE-3, requires an associated regulatory protein. *Proc. Natl. Acad. Sci. USA* 94: 3010-3015.
- Poulat, F., de Santa Barbara, P., Desclozeaux, M., Soullier, S., Moniot, B., Bonneaud, N., Boizet, B. and Berta, P. 1997. The human testis determining factor SRY bind a nuclear factor containing PDZ protein interaction domains. *J. Biol. Chem.* 272: 7167-7172.
- Reczek, D. and Bretscher, A. 2001. Identification of EPI64, a TBC/RabGAP domain-containing microvillar protein that binds to the first PDZ domain of EBP50 and E3KARP. *J. Cell Biol.* 153: 191-205.
- Itoh, T and Fukuda, M. 2006. Identification of EPI64 as a GTPase-activating protein specific for Rab27A. *J. Biol. Chem.* 281: 31823-31831.

CHROMOSOMAL LOCATION

Genetic locus: Tbc1d10a (mouse) mapping to 11 A1.

PRODUCT

EPI64 (m2): 293T Lysate represents a lysate of mouse EPI64 transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

APPLICATIONS

EPI64 (m2): 293T Lysate is suitable as a Western Blotting positive control for mouse reactive EPI64 antibodies. Recommended use: 10-20 µl per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

EPI64 (B-9): sc-376991 is recommended as a positive control antibody for Western Blot analysis of enhanced mouse EPI64 expression in EPI64 transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:
 1) Western Blotting: use m-IgGκ BP-HRP: sc-516102 or m-IgGκ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

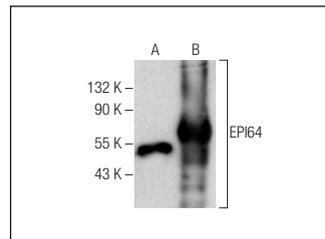

EPI64 (B-9): sc-376991. Western blot analysis of EPI64 expression in non-transfected: sc-117752 (A) and mouse EPI64 transfected: sc-125304 (B) 293T whole cell lysates.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.