

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



SOD-2 (h): 293T Lysate: sc-129809



The Power to Question

BACKGROUND

The superoxide dismutase family is composed of three metalloenzymes (SOD-1, SOD-2 and SOD-3) that catalyze the oxido-reduction of reactive oxygen species (ROS) such as superoxide anion. The SOD-2 precursor is a 222 amino acid protein that is encoded by nuclear chromatin, synthesized in the cytosol and imported post-translationally into the mitochondrial matrix. Unlike SOD-1, which is a homodimeric cytosolic Cu-Zn enzyme, SOD-2 is a homotetrameric manganese enzyme (also known as MnSOD) that functions in the mitochondrion. ROS are implicated in a wide range of degenerative processes, including Alzheimer's disease, Parkinson's disease and ischemic heart disease. Homozygous mutant mice, which lack SOD-2, exhibit dilated cardiomyopathy, accumulation of lipid in liver and skeletal muscle, metabolic acidosis, oxidative DNA damage and respiratory chain deficiencies in heart and skeletal muscle. Polymorphisms in the SOD-2 gene have also been implicated in nonfamilial, idiopathic, dilated cardiomyopathy in humans.

REFERENCES

- Wispe, J.R., et al. 1989. Synthesis and processing of the precursor for human mangano-superoxide dismutase. Biochem. Biophys. Acta 994: 30-36.
- Nishi, H., et al. 1995. DNA typing of HLA class II genes in Japanese patients with dilated cardiomyopathy. J. Mol. Cell. Cardiol. 27: 2385-2392.
- 3. Li, Y., et al. 1995. Dilated cardiomyopathy and neonatal lethality in mutant mice lacking manganese superoxide dismutase. Nat. Genet. 11: 376-381.
- Borgstahl, G.E., et al. 1996. Human mitochondrial manganese superoxide dismutase polymorphic variant Ile58Thr reduces activity by destabilizing the tetrameric interface. Biochemistry 35: 4287-4297.
- Hsieh, Y., et al. 1998. Probing the active site of human manganese superoxide dismutase: the role of glutamine 143. Biochemistry 37: 4731-4739.
- Melov, S., et al. 1998. A novel neurological phenotype in mice lacking mitochondrial manganese superoxide dismutase. Nat. Genet. 18: 159-163.
- 7. Melov, S., et al. 1999. Mitochondrial disease in superoxide disumtase 2 mutant mice. Proc. Natl. Acad. Sci. USA 96: 846-851.

CHROMOSOMAL LOCATION

Genetic locus: SOD2 (human) mapping to 6q25.3.

PRODUCT

SOD-2 (h): 293T Lysate represents a lysate of human SOD-2 transfected 293T cells and is provided as 100 μg protein in 200 μl SDS-PAGE buffer.

APPLICATIONS

SOD-2 (h): 293T Lysate is suitable as a Western Blotting positive control for human reactive SOD-2 antibodies. Recommended use: 10-20 μ l per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

SOD-2 (B-1): sc-133254 is recommended as a positive control antibody for Western Blot analysis of enhanced human SOD-2 expression in SOD-2 transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended: 1) Western Blotting: use m-lgG κ BP-HRP: sc-516102 or m-lgG κ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz MarkerTM Molecular Weight Standards: sc-2035, UltraCruz[®] Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

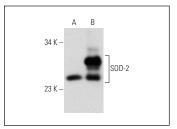

SOD-2 (B-1): sc-133254. Western blot analysis of SOD-2 expression in non-transfected: sc-117752 (A) and human SOD-2 transfected: sc-129809 (B) 293T whole cell lysates

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.

Santa Cruz Biotechnology, Inc. 1.800.457.3801 831.457.3801 Fax 831.457.3801 Europe +00800 4573 8000 49 6221 4503 0 www.scbt.com