



**SZABO
SCANDIC**

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic



METTL7B (h): 293T Lysate: sc-158726

BACKGROUND

METTL7B (methyltransferase-like protein 7B) is a 244 amino acid protein belonging to the methyltransferase superfamily. METTL7B is believed to have methyltransferase activity, wherein METTL7B catalyzes the transfer of a methyl group from one compound to another. The gene that encodes METTL7B maps to chromosome 12 which makes up about 4.5% of the human genome. A number of skeletal deformities are linked to chromosome 12 including hypochondrogenesis, achondrogenesis and Kniest dysplasia. Chromosome 12 is also home to a homeobox gene cluster which encodes crucial transcription factors for morphogenesis, and the natural killer complex gene cluster encoding C-type lectin proteins which mediate the NK cell response to MHC I interaction. Trisomy 12p leads to facial development defects, seizure disorders and a host of other symptoms varying in severity depending on the extent of mosaicism and is most severe in cases of complete trisomy.

REFERENCES

- Allen, T.L., Brothman, A.R., Carey, J.C. and Chance, P.F. 1996. Cytogenetic and molecular analysis in trisomy 12p. Am. J. Med. Genet. 63: 250-256.
- Yang, W. and Cole, W.G. 1998. Low basal transcripts of the COL2A1 collagen gene from lymphoblasts show alternative splicing of exon 12 in the Kniest form of spondyloepiphyseal dysplasia. Hum. Mutat. 1: S1-S2.
- Trowsdale, J., Barten, R., Haude, A., Stewart, C.A., Beck, S. et al. 2001. The genomic context of natural killer receptor extended gene families. Immunol. Rev. 181: 20-38.
- Zumkeller, W., Volleth, M., Muschke, P., Tönnies, H., Heller, A., Liehr, T., Wieacker, P. and Stumm, M. 2004. Genotype/phenotype analysis in a patient with pure and complete trisomy 12p. Am. J. Med. Genet. A 129: 261-264.
- Kelley, J., Walter, L. and Trowsdale, J. 2005. Comparative genomics of natural killer cell receptor gene clusters. PLoS Genet. 1: 129-139.
- Kemmer, L.A., Cowan, J.M., Hoffman, J.D. and Bianchi, D.W. 2006. The natural history of trisomy 12p. Am. J. Med. Genet. A 140: 695-703.

CHROMOSOMAL LOCATION

Genetic locus: METTL7B (human) mapping to 12q13.2.

PRODUCT

METTL7B (h): 293T Lysate represents a lysate of human METTL7B transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

APPLICATIONS

METTL7B (h): 293T Lysate is suitable as a Western Blotting positive control for human reactive METTL7B antibodies. Recommended use: 10-20 µl per lane.

Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.


METTL7B (D-2): sc-398626 is recommended as a positive control antibody for Western Blot analysis of enhanced human METTL7B expression in METTL7B transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:

- 1) Western Blotting: use m-IgG_κ BP-HRP: sc-516102 or m-IgG_κ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

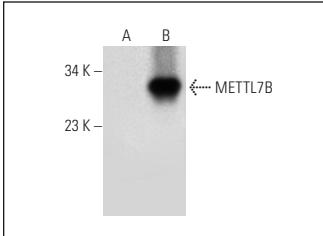

METTL7B (D-2): sc-398626. Western blot analysis of METTL7B expression in non-transfected: sc-117752 (**A**) and human METTL7B transfected: sc-158726 (**B**) 293T whole cell lysates.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.