



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

WASP (h): 293T Lysate: sc-170305

BACKGROUND

The Wiskott-Aldrich syndrome (WAS) is a disorder that results from a monogenic defect that has been mapped to the short arm of the X chromosome. WAS is characterized by thrombocytopenia, eczema, defects in cell-mediated and humoral immunity and a propensity for lymphoproliferative disease. The gene that is mutated in the syndrome encodes a proline-rich protein of unknown function designated WAS protein (WASP). A clue to WASP function came from the observation that T cells from affected males had an irregular cellular morphology and a disarrayed cytoskeleton suggesting the involvement of WASP in cytoskeletal organization. Close examination of the WASP sequence revealed a putative Cdc42/Rac interacting domain, homologous with those found in PAK65 and ACK. Subsequent investigation has shown WASP to be a true downstream effector of Cdc42.

REFERENCES

1. Remold-O'Donnell, E., et al. 1996. Defects in Wiskott-Aldrich syndrome blood cells. *Blood* 87: 2621-2631.
2. Stewart, D.M., et al. 1996. Studies of the expression of the Wiskott-Aldrich syndrome protein. *J. Clin. Invest.* 97: 2627-2634.
3. Symons, M., et al. 1996. Wiskott-Aldrich syndrome protein, a novel effector for the GTPase Cdc42Hs, is implicated in actin polymerization. *Cell* 84: 723-734.
4. Kolluri, R., et al. 1996. Direct interaction of the Wiskott-Aldrich syndrome protein with the GTPase Cdc42. *Proc. Natl. Acad. Sci. USA* 93: 5615-5618.
5. Schindelhauer, D., et al. 1996. Wiskott-Aldrich syndrome: no strict genotype-phenotype correlations but clustering of missense mutations in the amino-terminal part of the WASP gene product. *Hum. Genet.* 98: 68-76.

CHROMOSOMAL LOCATION

Genetic locus: WAS (human) mapping to Xp11.23.

PRODUCT

WASP (h): 293T Lysate represents a lysate of human WASP transfected 293T cells and is provided as 100 µg protein in 200 µl SDS-PAGE buffer.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

APPLICATIONS

WASP (h): 293T Lysate is suitable as a Western Blotting positive control for human reactive WASP antibodies. Recommended use: 10-20 µl per lane.


Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

WASP (B-9): sc-13139 is recommended as a positive control antibody for Western Blot analysis of enhanced human WASP expression in WASP transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).


RECOMMENDED SUPPORT REAGENTS

To ensure optimal results, the following support reagents are recommended:
1) Western Blotting: use m-IgGκ BP-HRP: sc-516102 or m-IgGκ BP-HRP (Cruz Marker): sc-516102-CM (dilution range: 1:1000-1:10000), Cruz Marker™ Molecular Weight Standards: sc-2035, UltraCruz® Blocking Reagent: sc-516214 and Western Blotting Luminol Reagent: sc-2048.

DATA

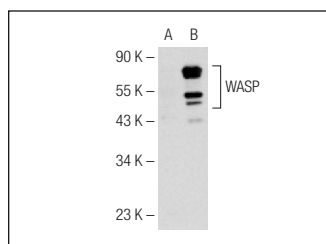

WASP (B-9): sc-13139. Western blot analysis of WASP expression in non-transfected: sc-117752 (A) and human WASP transfected: sc-170305 (B) 293T whole cell lysates.

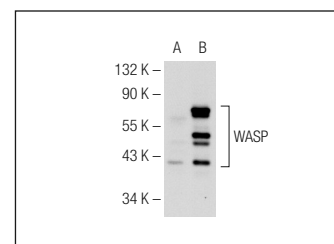

WASP (D-1): sc-5300. Western blot analysis of WASP expression in non-transfected: sc-117752 (A) and human WASP transfected: sc-170305 (B) 293T whole cell lysates.

RESEARCH USE

For research use only, not for use in diagnostic procedures.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.