



**SZABO
SCANDIC**

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic



AChR β 1 (h): 293T Lysate: sc-171924

BACKGROUND

Members of the ligand-gated ion channel receptor family are characterized by their fast transmitting response to neurotransmitters. Two important members of this family are the nicotinic acetylcholine and glutamate receptors, both of which are composed of five homologous subunits forming a transmembrane aqueous pore. These transmembrane receptors change conformation in response to their cognate neurotransmitter. Nicotinic acetylcholine receptors (AChRs) are found at the postsynaptic membrane of the neuromuscular junction and bind acetylcholine molecules, allowing ions to move through the pore. Glutamate receptors are found in the postsynaptic membrane of cells in the central nervous system. The activity that is generated at the synapse by the binding of acetylcholine is terminated by acetylcholinesterase, an enzyme that rapidly hydrolyzes acetylcholine.

REFERENCES

- Alkondon, M., Rao, K.S. and Albuquerque, E.X. 1988. Acetylcholinesterase reactivators modify the functional properties of the nicotinic acetylcholine receptor ion channel. *J. Pharmacol. Exp. Ther.* 245: 543-556.
- Betz, H. 1990. Ligand-gated ion channels in the brain: the amino acid receptor superfamily. *Neuron* 5: 383-392.
- Baenziger, J.E., Miller, K.W., McCarthy, M.P. and Rothschild, K.J. 1992. Probing conformational changes in the nicotinic acetylcholine receptor by Fourier transform infrared difference spectroscopy. *Biophys. J.* 62: 64-66.
- Daw, N.W., Stein, P.S. and Fox, K. 1993. The role of NMDA receptors in information processing. *Annu. Rev. Neurosci.* 16: 207-222.
- Unwin, N. 1993. Neurotransmitter action: opening of ligand-gated ion channels. *Cell* 72: 31-41.
- Stevens, C.F. 1993. Quantal release of neurotransmitter and long-term potentiation. *Cell* 72: 55-63.
- Sargent, P.B. 1993. The diversity of neuronal nicotinic acetylcholine receptors. *Annu. Rev. Neurosci.* 16: 403-443.
- Ramirez-Latorre, J., Yu, C.R., Qu, X., Perin, F., Karlin, A. and Role, L. 1996. Functional contributions of α 5 subunit to neuronal acetylcholine receptor channels. *Nature* 380: 347-351.

CHROMOSOMAL LOCATION

Genetic locus: CHRN1 (human) mapping to 17p13.1.

PRODUCT

AChR β 1 (h): 293T Lysate represents a lysate of human AChR β 1 transfected 293T cells and is provided as 100 μ g protein in 200 μ l SDS-PAGE buffer.

STORAGE

Store at -20° C. Repeated freezing and thawing should be minimized. Sample vial should be boiled once prior to use. Non-hazardous. No MSDS required.

RESEARCH USE

For research use only, not for use in diagnostic procedures.


APPLICATIONS

AChR β 1 (h): 293T Lysate is suitable as a Western Blotting positive control for human reactive AChR β 1 antibodies. Recommended use: 10-20 μ l per lane.

Control 293T Lysate: sc-117752 is available as a Western Blotting negative control lysate derived from non-transfected 293T cells.

AChR β 1 (118): sc-65787 is recommended as a positive control antibody for Western Blot analysis of enhanced human AChR β 1 expression in AChR β 1 transfected 293T cells (starting dilution 1:100, dilution range 1:100-1:1,000).

DATA

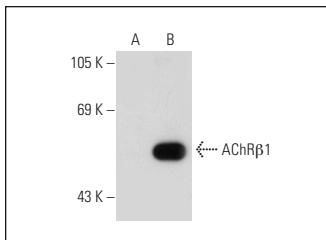

AChR β 1 (118): sc-65787. Western blot analysis of AChR β 1 expression in non-transfected: sc-117752 (**A**) and human AChR β 1 transfected: sc-171924 (**B**) 293T whole cell lysates.

PROTOCOLS

See our web site at www.scbt.com for detailed protocols and support products.