

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten! See the following pages for more information!



## Lieferung & Zahlungsart

siehe unsere Liefer- und Versandbedingungen

## Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

linkedin.com/company/szaboscandic in



# Vmn2r114 HDR Plasmid (m): sc-437088-HDR



The Power to Question

#### **BACKGROUND**

DNA containing double-strand breaks (DSB) created by the CRISPR/Cas9 system can be repaired by either the non-homologous end-joining (NHEJ) or the homology-directed repair (HDR) pathway (1,2,3). The NHEJ repair pathway introduces non-specific insertions or deletions at the cleavage site, whereas the HDR pathway allows for precise gene editing at the DSB site (1,2,3). Target-specific HDR Plasmids provide a DNA repair template for a DSB and, when co-transfected with CRISPR/Cas9 KO Plasmids, enable the insertion of specific selection markers where Cas9-induced DNA cleavage has occurred (1,2). The HDR plasmid can incorporate a Red Fluorescent Protein (RFP) gene to visually confirm transfection and an antibiotic resistance gene (puromycin) for selection of cells containing a successful CRISPR/Cas9 double-strand break. The puromycin resistance and RFP encoding genes are flanked by two LoxP sites that are recognized by the Cre Vector, which can be used to later remove these selection markers from the genomic DNA (4,5).

#### **REFERENCES**

- Mali, P., et al. 2013. RNA-guided human genome engineering via Cas9. Science 339: 823-826.
- Ran, F.A., et al. 2013. Genome engineering using the CRISPR-Cas9 system. Nat. Protoc. 8: 2281-2308.
- 3. Hsu, P., et al. 2014. Development and applications of CRISPR-Cas9 for genome editing. Cell 157: 1262-1278.
- Ma, Y. 2014. Generation of eGFP and Cre knockin rats by CRISPR/Cas9. FEBS J. 281: 3779-3790.
- 5. Ma, Y., et al. 2014. Generating rats with conditional alleles using CRISPR/Cas9. Cell Res. 24: 122-125.

#### **CHROMOSOMAL LOCATION**

Genetic locus: Vmn2r114 (mouse) mapping to 17 A3.3.

#### **PRODUCT**

Vmn2r114 HDR Plasmid (m) consists of a pool of 2-3 plasmids, each containing a homology-directed DNA repair (HDR) template corresponding to the cut sites generated by the Vmn2r114 CRISPR/Cas9 KO Plasmid (m): sc-437088. Each HDR template contains two 800 bp homology arms designed to specifically bind to the genomic DNA surrounding the corresponding Cas9-induced double-strand DNA break site. Each vial contains 20 µg of lyophilized HDR Plasmid DNA. Suitable for up to 20 transfections.


#### STORAGE AND RESUSPENSION

Store lyophilized plasmid DNA at  $4^{\circ}$  C with desiccant. Stable for at least one year from the date of shipment. Once resuspended, store at  $4^{\circ}$  C for short term storage or  $-20^{\circ}$  C for long-term storage. Avoid repeated freeze thaw cycles.

Resuspend lyophilized plasmid DNA in 200  $\mu$ l of the provided ultrapure, sterile, DNase-free water. Resuspension of the plasmid DNA makes a 0.1  $\mu$ g/ $\mu$ l solution in a 10 mM TRIS EDTA, 1 mM EDTA buffered solution.

#### **APPLICATIONS**

Vmn2r114 HDR Plasmid (m) is recommended for co-transfection with Vmn2r114 CRISPR/Cas9 KO Plasmid (m): sc-437088 and designed for repair of the site-specific Cas9-induced DNA cleavage within the Vmn2r114 (mouse) gene. During repair, the Vmn2r114 HDR Plasmid (m) incorporates a puromycin resistance gene to enable selection of stable knockout (KO) cells and an RFP gene to visually confirm transfection.

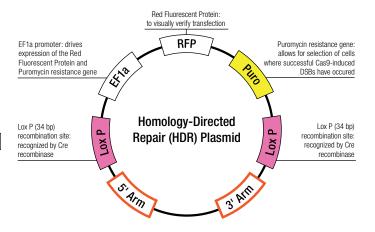

#### **SUPPORT REAGENTS**

For optimal reaction efficiency with HDR Plasmids, Santa Cruz Biotechnology's UltraCruz® Transfection Reagent: sc-395739 (0.2 ml), Plasmid Transfection Medium: sc-108062 (20 ml) and L-755,507: sc-204045 (10 mg) are recommended. Puromycin dihydrochloride: sc-108071 (25 mg) is recommended for selection. Cre Vector: sc-418923 (20  $\mu g$  in 20  $\mu l$ ) is also available for the optional removal of the puromycin resistance gene inserted during homology-directed repair.

#### **RESEARCH USE**

The purchase of this product conveys to the buyer the nontransferable right to use the purchased amount of the product and all replicates and derivatives for research purposes conducted by the buyer in his laboratory only (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party, or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes.

#### **PROTOCOLS**

See our web site at www.scbt.com or our catalog for detailed protocols and support products.

Santa Cruz Biotechnology, Inc. 1.800.457.3801 831.457.3800 fax 831.457.3801 Europe +00800 4573 8000 49 6221 4503 0 www.scbt.com