



# SZABO SCANDIC

Part of Europa Biosite

## Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!  
See the following pages for more information!



### Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

### Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

### SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

[mail@szabo-scandic.com](mailto:mail@szabo-scandic.com)

[www.szabo-scandic.com](http://www.szabo-scandic.com)

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

# mouse anti-goat IgG-CFL 680: sc-516245

## BACKGROUND

Santa Cruz Biotechnology's high quality, well characterized monoclonal secondary antibodies are available conjugated to either an enzyme, biotin or fluorophore for use in a variety of antibody-based applications, including Western blotting, immunostaining and flow cytometry. Santa Cruz secondary antibodies are commonly affinity purified against immobilized whole IgG isotypes, including IgG<sub>1</sub>, IgG<sub>2a</sub>, IgG<sub>2b</sub>, IgG<sub>3</sub> and IgG<sub>4</sub>. Monoclonal secondary antibodies are available conjugated to HRP for Western blotting (WB) and immunohistochemistry (IHC); (CM) or Cruz Marker form of HRP conjugated secondary antibodies are suitable for use with our Cruz Marker™ molecular weight standards; FITC (fluorescein isothiocyanate), PE (phycoerythrin), R (TRITC: tetramethyl rhodamine isothiocyanate), TR (Texas Red®), PerCP (peridinin chlorophyll protein complex), PerCP-Cy5.5 (peridinin chlorophyll protein complex with cyanin-5.5), and CruzFluor™ (488, 555, 594 and 647) for immunofluorescence (IF), immunohistochemistry (IHC) and flow cytometry (FCM); B (biotin) for immunohistochemistry (IHC); AP (alkaline phosphatase) for Western blotting (WB); and CruzFluor® 680 and 790 for near-infrared (NIR) Western blotting (WB), immunofluorescence (IF), immunohistochemistry (IHC) and flow cytometry (FCM).

## SOURCE

mouse anti-goat IgG-CFL 680 is a pre-adsorbed, affinity purified secondary antibody raised in mouse against goat IgG and conjugated to CruzFluor™ 680 (CFL 680).

## PRODUCT

Each vial contains 200 µg mouse IgG (pre-adsorbed with mouse and human IgG) in 0.5 ml of PBS containing 0.1% gelatin and 0.1% sodium azide.

## APPLICATIONS

mouse anti-goat IgG-CFL 680 is recommended for detection of goat IgG by NIR Western blotting (starting dilution: 1:100, dilution range: 1:100-1:400), immunofluorescence staining (starting dilution: 1:100, dilution range: 1:100-1:400), immunohistochemical staining (starting dilution: 1:100, dilution range: 1:100-1:400) and flow cytometry (0.5-1 µg per 1 x 10<sup>6</sup> cells). Optimal dilution to be determined by titration.

## STORAGE

Store at 4° C, **\*\*DO NOT FREEZE\*\***. Stable for one year from the date of shipment. Non-hazardous. No MSDS required.

## PROTOCOLS

See our web site at [www.scbt.com](http://www.scbt.com) for detailed protocols and support products.

## RECOMMENDED SUPPORT PRODUCTS

### A. TISSUE CULTURE CELLS

- CrystalCruz™ Cover Glasses, 22 x 50 mm, precleaned: sc-24975
- PBS (Phosphate Buffered Saline), powder, 1 packet: sc-24947
- Hydrogen Peroxide, 30% solution, 100 ml: sc-203336


### B. FROZEN TISSUE SECTIONS

- Organo/Limonene Mount, non-toxic alternative to Permount, 100 ml: sc-45087
- UltraCruz® Mounting Medium, aqueous-based, 10 ml: sc-24941
- ImmunoHistoMount, aqueous-based mounting medium, 30 ml: sc-45086
- Immuno In Situ Mount, for use with *in situ* hybridization, 30 ml: sc-45088


### C. FORMALIN-FIXED, PARAFFIN-EMBEDDED TISSUE SECTIONS

- Paraffin, for the preparation of tissue samples for staining, 500 g: sc-286633
- Xylenes, mixed isomers with ethylbenzene, 500 ml: sc-237422
- Hematoxylin, Gill's Formulation #2; nuclear counter stain, 100 ml: sc-24973


## DATA

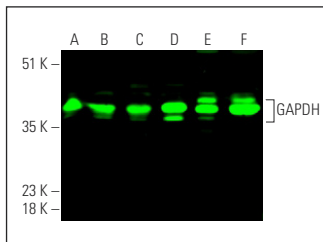


GAPDH (V-18): sc-20357. Near-infrared western blot analysis of GAPDH expression in HeLa (A), Hep G2 (B), A549 (C), JAR (D), NIH/3T3 (E) and KNRK (F) whole cell lysates. Detection reagent used: mouse anti-goat IgG-CFL 680: sc-516245.

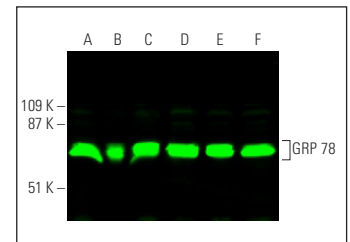

GRP 78 (N-20): sc-1050. Near-infrared western blot analysis of GRP 78 expression in Jurkat (A), HeLa (B), Hep G2 (C), JAR (D), NIH/3T3 (E) and KNRK (F) whole cell lysates. Detection reagent used: mouse anti-goat IgG-CFL 680: sc-516245.

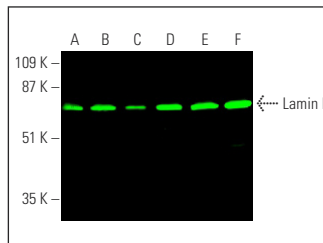

Lamin B (M-20): sc-6217. Near-infrared western blot analysis of Lamin B expression in Jurkat (A), Y79 (B), HeLa (C), MOLT-4 (D), CCRF-CEM (E) and F9 (F) whole cell lysates. Detection reagent used: mouse anti-goat IgG-CFL 680: sc-516245.

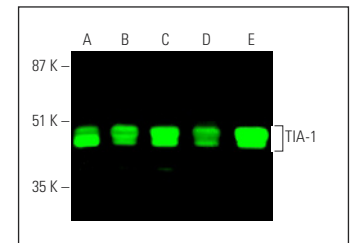

TIA-1 (C-20): sc-1751. Near-infrared western blot analysis of TIA-1 expression in Jurkat (A), BJAB (B), K-562 (C), MOLT-4 (D) and HuT 78 (E) whole cell lysates. Detection reagent used: mouse anti-goat IgG-CFL 680: sc-516245.

## RESEARCH USE

For research use only, not for use in diagnostic procedures.