



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7

mail@szabo-scandic.com

www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 


Cresyl Violet Acetate (1.0%, aqueous)

Description: Cresyl Violet Acetate (1.0%, aqueous) is designed for staining Nissl substance in neurons on formalin fixed, paraffin-embedded tissue. Nissl granules are purple, nuclei of neuroglia and endothelial cells are slightly bluer than Nissl granules.

Nissl Substance: Violet
Neuroglia Nuclei: Violet to Dark Blue

Uses/Limitations: Not to be taken internally.
For Research use only.
Histological applications.
Do not use if reagent become cloudy.
Do not use past expiration date.
Use caution when handling reagents.
Non-Sterile.

Control Tissue: Cerebral Cortex

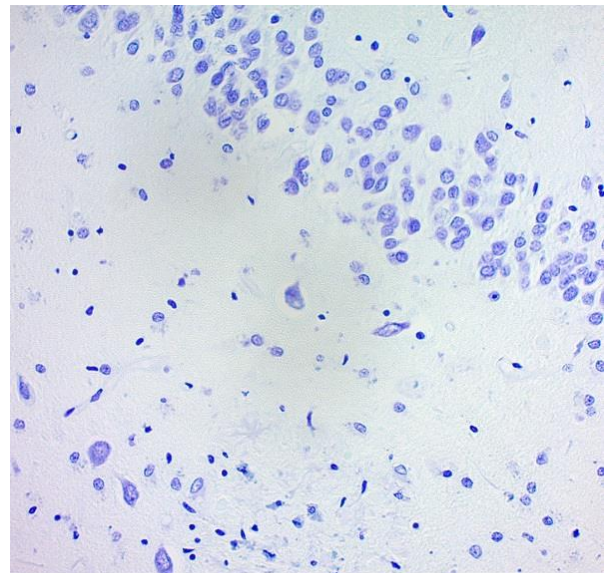

Cresyl Violet Acetate (1.0%, aqueous) diluted 1:1 with Deionized water staining of human brain at 40X magnification

Availability/Contents:

<u>Item #</u>	<u>Volume</u>	<u>Storage</u>
CEB500	500 ml	2-8° C
CEB999	1000 ml	2-8° C

Precautions: Avoid contact with skin and eyes.
Follow all Federal, State, and local regulations regarding disposal.

Procedure:

Working Solutions: If preferred, reagent may be diluted up to 1:10 with Deionized water just before use.

1. Deparaffinize sections and hydrate to distilled water.
2. Apply Cresyl Violet Acetate solution (or working solution) to tissue for 3-5 minutes
3. Rinse quickly in 1 change of distilled water.

Storage: 2° C



8° C



ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.

**C
P**

Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands

4. Dehydrate quickly in absolute alcohol. (Alcohol may remove stain from tissue over time)
5. Clear in Xylene or Xylene substitute and mount in synthetic resin.

References:

1. Sheenan, D.C., Hrapchak, B.B. Theory and Practice of Histotechnology, 2nd Edition. Battelle Press, Columbus, OH. Page 262-264. 1980
2. Kluver, H., Barrera, E.A. A Method for the combined staining of cells and fibers in the nervous system. Journal of Neuropathology and Experimental Neurology, 1953, 12: pages 400-403.
3. Margaret M. Powers & George Clark (1955) An Evaluation of Cresyl Echt Violet Acetate as a Nissl Stain, Stain Technology, 30:2, 83-88,

Storage: 2° C



8° C

ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.**C
P**Emergo Europe
Prinsessegracht 20
2514 AP The Hague, The Netherlands