



SZABO SCANDIC

Part of Europa Biosite

Produktinformation



Forschungsprodukte & Biochemikalien



Zellkultur & Verbrauchsmaterial



Diagnostik & molekulare Diagnostik



Laborgeräte & Service

Weitere Information auf den folgenden Seiten!
See the following pages for more information!



Lieferung & Zahlungsart

siehe unsere [Liefer- und Versandbedingungen](#)

Zuschläge

- Mindermengenzuschlag
- Trockeneiszuschlag
- Gefahrgutzuschlag
- Expressversand

SZABO-SCANDIC HandelsgmbH

Quellenstraße 110, A-1100 Wien

T. +43(0)1 489 3961-0

F. +43(0)1 489 3961-7


mail@szabo-scandic.com


www.szabo-scandic.com

[linkedin.com/company/szaboscandic](https://www.linkedin.com/company/szaboscandic) 

CD31 / PECAM-1 (Endothelial Cell Marker); Clone 1A10 (Concentrate)


Availability/Contents:	<u>Item #</u>	<u>Volume</u>
	RA0260-C.5	0.5 ml
Description:		
Species:	Mouse	
Immunogen:	Prokaryotic recombinant fusion protein corresponding to a portion of the extracellular domain downstream of the signal sequence of the CD31 (PECAM-1).	
Clone:	1A10	
Isotype:	IgG1, kappa	
Entrez Gene ID:	5175 (Human); 18613 (Mouse)	
Hu Chromosome Loc.:	17q23.3	
Synonyms:	EndoCAM; PECA1; Platelet Endothelial Cell Adhesion Molecule 1; GPIIA'	
Mol. Weight of Antigen:	~100kDa (endothelium) and ~130kDa (platelets)	
Format:	Tissue culture supernatant with 0.05% Azide.	
Specificity:	Anti-CD31 has shown to be highly specific and sensitive for vascular endothelial cells. Staining of nonvascular tumors (excluding hematopoietic neoplasms) is rare. Anti-CD31 reacts with normal, benign, and malignant endothelial cells which make up blood vessel lining.	
Background:	CD31 (PECAM-1) is a transmembrane glycoprotein member of the immunoglobulin supergene family of adhesion molecules. CD31 is expressed by stem cells of the hematopoietic system and is primarily used to identify and concentrate these cells for experimental studies as well as for bone marrow transplantation. The level of CD31 expression can help to determine the degree of tumor angiogenesis, and a high level of CD31 expression may imply a rapidly growing tumor and potentially be a predictor of tumor recurrence.	
Species Reactivity:	Human. Others not known.	
Positive Control:	Tonsil, Angiosarcoma.	
Cellular Localization:	Cell surface and cytoplasmic	
Titer/ Working Dilution:	Immunohistochemistry (Frozen and Formalin-fixed): 1:50-1:100	
	Flow Cytometry:	5-10 µl/million cells
	Immunofluorescence:	1:50-1:100
	Western Blotting:	1:100-1:200
	Immunoprecipitation:	5-10 µl/500µg protein lysate
Microbiological State:	This product is not sterile.	

Storage: 2° C  8° C

 ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.

  EmergoEurope (31)(0) 70 345-8570
Molsnstraat 15
2513 BH Hague, The Netherlands

Uses/Limitations: Not to be taken internally.
For Research Use Only.
This product is intended for qualitative immunohistochemistry with normal and neoplastic formalin-fixed, paraffin-embedded tissue sections, to be viewed by light microscopy.
Do not use if reagent becomes cloudy.
Do not use past expiration date.
Non-Sterile.

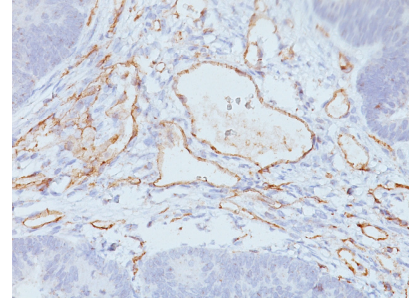

Ordering Information and Current Pricing at www.scytek.com

Formalin-fixed, paraffin-embedded tonsil stained with CD31; Clone 1A10.

Procedure:


- Tissue Section Pretreatment (Required):** Staining of formalin fixed, paraffin embedded tissue sections is significantly enhanced by pretreatment with Citrate Plus (ScyTek catalog# CPL500).
- Primary Antibody Incubation Time:** We suggest an incubation period of 30 minutes at room temperature. However, depending upon the fixation conditions and the staining system employed, optimal incubation should be determined by the user.
- Visualization:** For maximum staining intensity we recommend the “UltraTek HRP Anti-Polyvalent Lab Pack” (ScyTek catalog# UHP125, see IFU for instructions) combined with the “DAB Chromogen/Substrate Bulk Pack (High Contrast)” (ScyTek catalog# ACV500, see IFU for instructions).


Precautions: Contains Sodium Azide as a preservative (0.09% w/v).
Do not pipette by mouth.
Avoid contact of reagents and specimens with skin and mucous membranes.
Avoid microbial contamination of reagents or increased nonspecific staining may occur.
This product contains no hazardous material at a reportable concentration according to U.S. 29 CFR 1910.1200, OSHA Hazardous Communication Standard and EC Directive 91/155/EC.

References:

- Stockinger H, Gadd S J, Eher R, et al. Molecular characterization and functional analysis of the leukocyte surface protein CD31. *Journal of Immunology*. 145 (11): 3889-3897 (1990).
- Stockinger H, Schreiber W, Majdic O, et al. Phenotype of human T cells expressing CD31, a molecule of the immunoglobulin supergene family. *Immunology*. 75 (1): 53-58 (1992).


Warranty: No products or “Instructions For Use (IFU)” are to be construed as a recommendation for use in violation of any patents. We make no representations, warranties or assurances as to the accuracy or completeness of information provided on our IFU or website. Our warranty is limited to the actual price paid for the product. ScyTek Laboratories, Inc. is not liable for any property damage, personal injury, time or effort or economic loss caused by our products. Immunohistochemistry is a complex technique involving both histological and immunological detection methods. Tissue processing and handling prior to immunostaining can cause inconsistent results. Variations in fixation and embedding or the inherent nature of the tissue specimen may cause variations in results. Endogenous peroxidase activity or pseudoperoxidase activity in erythrocytes and endogenous biotin may cause non-specific staining depending on detection system used.

Storage: 2° C  8° C



ScyTek Laboratories, Inc.
205 South 600 West
Logan, UT 84321
U.S.A.

CE

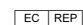 EmergoEurope (31)(0) 70 345-8570
Molsnstraat 15
2513 BH Hague, The Netherlands